VERMEIREN

V300 - 30°

INSTRUCTION MANUAL
MODE D'EMPLOI
GEBRUIKSAANWIJZING
GEBRAUCHSANWEISUNG
ISTRUZIONI PER L'USO
MANUAL DE INSTRUCCIONES
INSTRUKCJA OBSŁUGI

V300 30° - Additions to manual V300

1 Product description (refer to manual V300 paragraph 1)

1.1 Technical specifications (add to manual V300 paragraph 1.1)


Technical terms below are valid for the wheelchair in standard settings. If other footrest / arm support or other accessories are used, the tabulated values will change.

Support of other accessories are used, the tabulated values will change.

Make	Vermeiren					
Address	Vermeirenplein 1/15, B-2920 Kalmthout					
Type	Manual wheelchair					
Model	V300 30°					
Maximum occupant mass	130 kg					
Description	Min.			Max.		
Overall length with legrest	1040 mm					
Effective seat width	390 mm	420 mm	440 mm	460 mm	480 mm	500 mm
Overall width (depends on the seat width)	590 mm	620 mm	640 mm	660 mm	680 mm	700 mm
Folded length	1020 mm					
Folded width	310 mm					
Folded height	910 mm					
Total mass	17,7 kg					
Mass of heaviest part	10 kg					
Static stability downhill	10° (in standard configuration)					
Static stability uphill	2°					
Static stability sideways	18° (in standard configuration)					
Obstacle climbing	60 mm					
Seat plane angle	0°			10°		
Effective seat depth	440 mm			520 mm		
Seat surface height at front edge	440 mm			530 mm		
Backrest angle	5°			35°		
Backrest height	400 mm					
Distance between footrest and seat	420 mm					
Angle between seat and footrest	5°			25°		
Distance between armrest and seat	220 mm			240 mm		
Front location of armrest structure	410 mm					
Handrim diameter	535 mm					
Horizontal location of axle (deflection)	1,27 mm			61,04 mm		
Minimum turning radius	1500 mm					
Diameter Krypton PU Rear wheels	22"			24"		
Tyre pressure, rear (driving) wheels	Max. 3.5 bar					
Diameter Krypton PU steering wheels	200 mm					
Tyre pressure, steering wheels	Max. 2.5 bar					
Storage and use temperature	+ 5 °C			+ 41 °C		
Storage and use humidity	30%			70%		
We reserve the right to introduce technical changes. Measurement tolerance ± 15 mm / 1,5 kg / °						

Table 1: Technical specifications V300 30°

1.2 Drawing (change paragraph 1.3 manual V300)

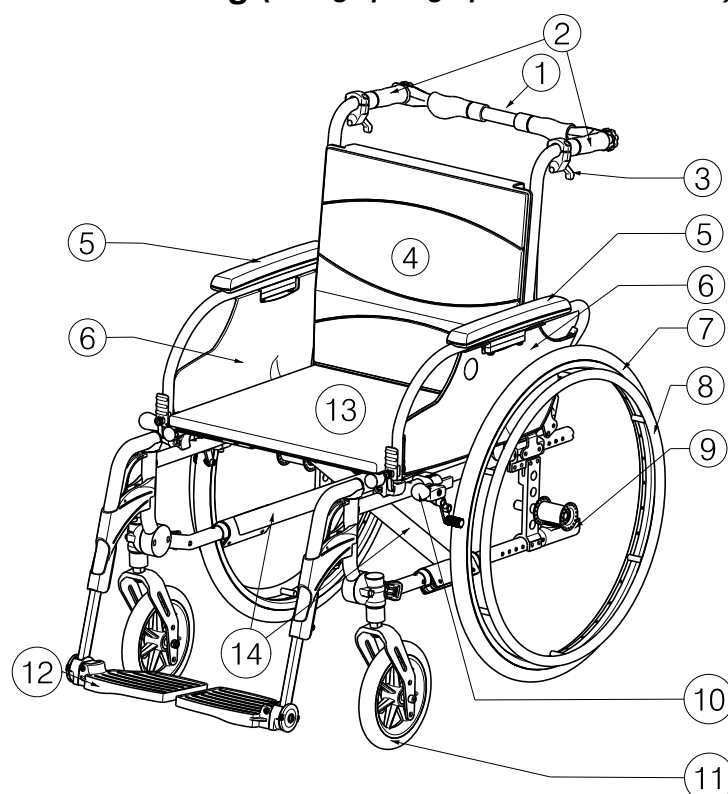- 1 = Push bar
- 2 = Handles
- 3 = Lever back inclination
- 4 = Backrest
- 5 = Armrests
- 6 = Arm supports
- 7 = Driving wheels (rear wheels)
- 8 = Hand rims
- 9 = Tip cap
- 10 = Brakes
- 11 = Steering wheels (front wheels)
- 12 = Footrests
- 13 = Seat
- 14 = Cross

1.3 Accessories (add to manual V300 paragraph 1.4)

The following accessories are available for the V300 30°:

- Drum brakes (B74)

1.4 Explanation of symbols (add to manual V300 paragraph 1.6)


Not crash-tested

2 Use (refer to manual V300 paragraph 2)


2.1 Backrest inclination (add to manual V300)

- ⚠ **WARNING:** Risk of tipping over – Be aware that the stability decreases when backrest is adjusted backwards.
- ⚠ **WARNING:** Risk of injury – Apply the parking brakes before adjusting the backrest.
- ⚠ **CAUTION:** Chance of pinching – Do not place your fingers, clothes between the adjustment mechanism.

With the model V300 30° it is possible to set the backrest backwards with a maximum inclination of 30° (4 different positions: steps 10°). Be sure that the patient sits in the chair when the attendant shall perform the back adjusting and that the wheelchair does not tip over.

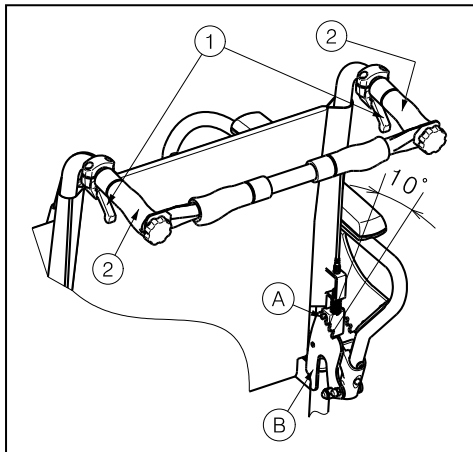
1. Pull handle ① on both sides towards the grab handle ② to adjust the backrest.
2. The pen ① comes out the groove of the toothed plate ②.
3. Pull both back tubes backwards continuously to the desired position (4 different positions: 5° - 15° - 25° - 35°).
4. Check that the pen ① is well fixed in the groove of the toothed plate ②.

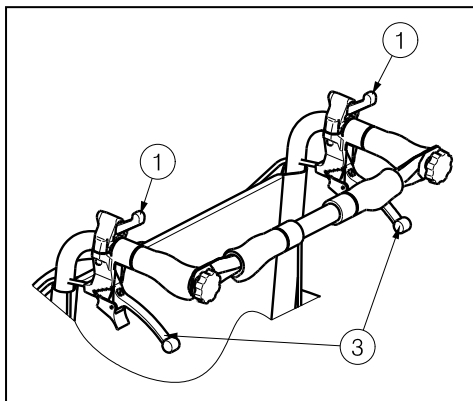If your wheelchair features drum brakes, the handles ③ are used as brakes.

2.2 Mounting or removing of push bar *(add to manual V300)*

- ⚠ **WARNING:** Risk of injury - Make sure that buttons ① is tightened.
- ⚠ **WARNING:** Risk of injury - When equipped with the push bar, the push bar shall always be mounted (→ to increase the stability of the back).

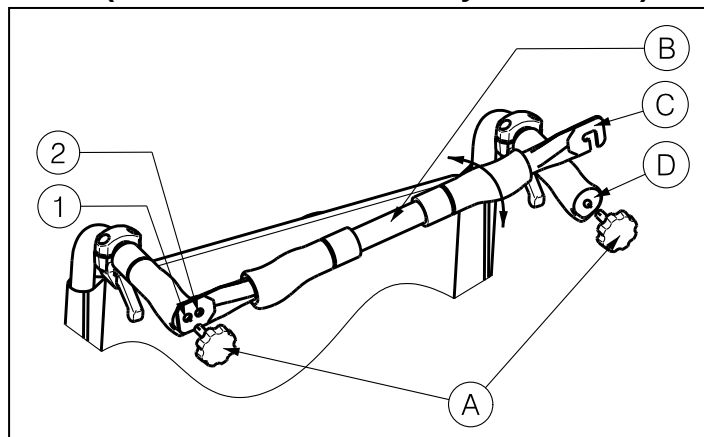The mounting of the push bar ② is done as follows:

1. Mount the push bar ② with the starknob ① in the left hand-grip.
2. Swing the push bar ② to the right.
3. Mount hook ③ over the right hand-grip ④ as shown on the drawing. Use holes ① or ② to adjust the tension of the backrest.

To take off the Push bar:

1. Loosen the starknob ① on the right hand-grip.
2. Remove hook ③ from the right hand-grip ④.

2.3 **Transport in the car** (add to manual V300 paragraph 2.14)

⚠ DANGER: Risk of injury – The wheelchair is not suited for use as a seat in a motor vehicle.

Never use your wheelchair as a seat in an automobile or other vehicle. The wheelchair shall be marked with following symbol.


3 **Installation and adjustment** (refer to manual V300 paragraph 3)

3.1 **Manner of delivery** (add to manual V300 paragraph 3.2)

The Vermeiren V300 30° shall be delivered with:

- 1 frame with arm supports, rear and front wheels (standard delivery seat height: 500 mm, seat angle: 5°)
- 1 pair footrests
- Back adjusting mechanism
- Tools
- Manual
- Accessories

3.2 **Adjusting the stability and the maneuverability** (add to manual V300 paragraph 3.8)

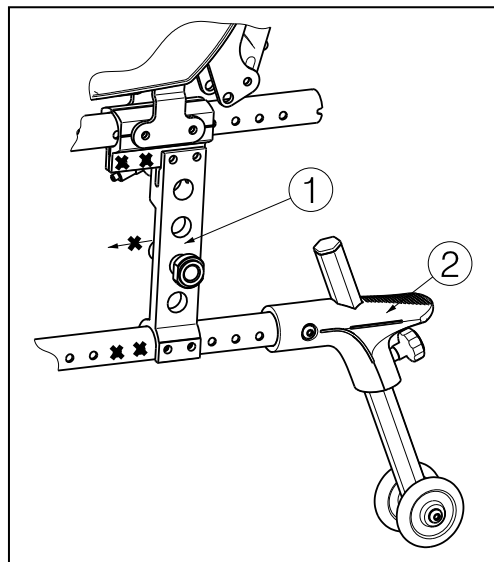Axle block ① in most backwards position, anti-tipping ② is standard available.

4 **Maintenance** (refer to manual V300 paragraph 4)

4.1 **Regular Maintenance** (add to manual V300 paragraph 4.1)

- Every 8 weeks: inspection and possibly lubricating or adjusting from
 - Backrest inclination system

4.2 **Inspection** (add to manual V300 paragraph 4.4)

- Check the functioning of the backrest inclination (locking, load, deformation, wear caused by loads).

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Vermeiren Nederland B.V.